

yuwell

Classical heritage of 1 million units

Homecare oxygen concentrator YU300



MEDICAL
BACKGROUND

Professional
medical company

GLOBAL
LEADER

Sales volume has been
1.1 million units

OXYGEN
GENERATOR

Supplying oxygen generators
for 130 hospitals


LISTED
COMPANY

Company's assets
reached 321 million USD

If you also got these feelings..



Insomnia, yawning,
head pain, cold feet

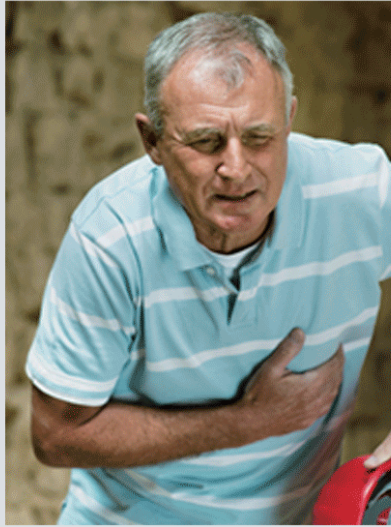

Bad breath, fast heartbeat



Pregnant woman who suffering
from breathing problem




Nervous, sleepy, tired

Attention




Medical research shows that Hypoxia may cause breathing problem, bad breath, dry skin, sleeping problem, nervous, high blood pressures, immun system drop down, etc.

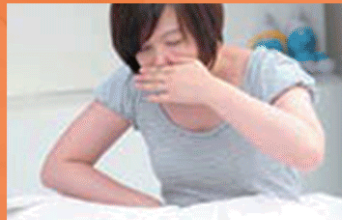

Serious hypoxia may cause more other diseases.
It is not in the alarmist

Who needs oxygen therapy?

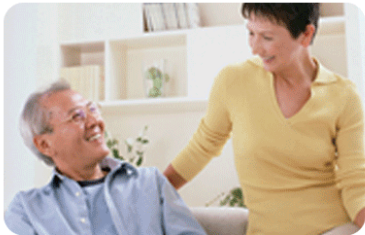

Aged over 60



Pregnant woman



Office workers,
white-collar



Students

Oxygen is good for cardiovascular, respiratory diseases, and is suitable for tracheitis, bronchitis, asthma, lung disease, respiratory failure, cerebral thrombosis, cerebral infarction, cerebral atrophy or hypoxia-induced headaches, insomnia or drowsiness, coronary heart disease, angina pectoris, myocardial infarction, hypertension, atherosclerosis, heart failure and other rehabilitation.

How to oxygenating?

Method 1

Keeping healthy life-style persevering.



1. Insist on exercising.



2. Rest regularly , getting plenty of sleep.



3. Being closer to nature.

Method 2

In one step

20 minutes for oxygen therapy everyday!



Oxygen **BAR** in your home

Bring you a new life-style with oxygen.

YU300
HOME CARE OXYGEN



Portable

Only 8kg, easy to carry, supplied with vehicle power, can be used in a car!



The pursuit of excellence is our
ultimate goal!

YU300 Homecare oxygen concentrator,
the essence of 1.1 million units.
Level up in technology, design and
appearance.


1.

Stable oxygen concentration

- 1 Stable: continuous working for 48 hours with no failure
Other homecare oxygen concentrator can not continuous work for more than 12 hours




Oxygen concentration

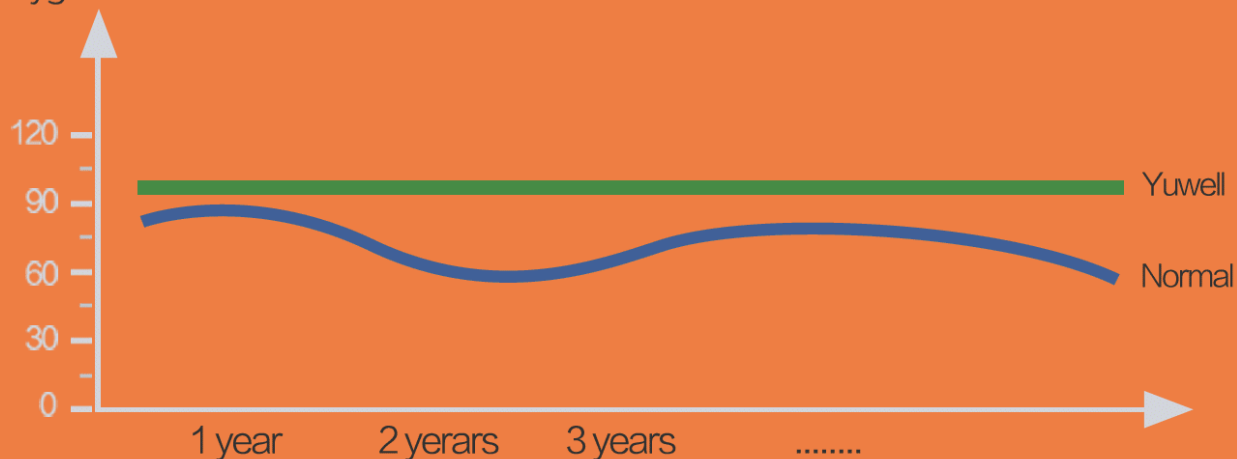

- 2 Oxygen concentration: 0 ~ 1L/min: $\geq 90\%$



Oxygen concentration

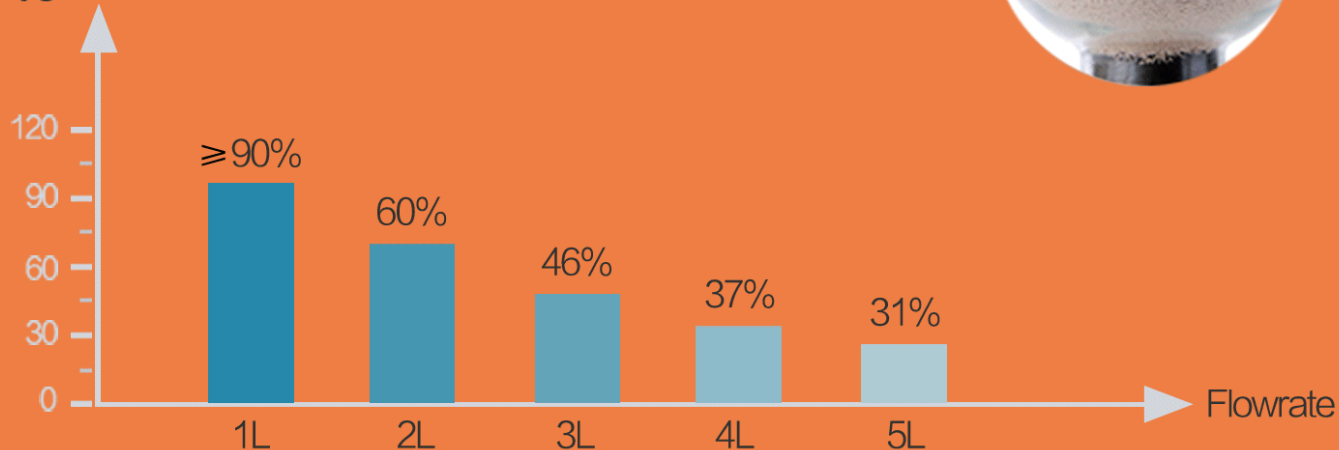

2.

Ultra silent

Noise $\leq 40\text{dB(A)}$



dB

0–20

Very quiet

dB

20–35

Quiet

≤ 40 dB



dB

40–50

Normal

dB

50–60

Noisy

3. Easy Shuttle



4.

Intelligent surveillance



Failure alarm



Filter cleaning
alarm



Real-time
oxygen concentration



5.


Infrared Remote Control System



6.

High reverse display LCD screen

The display screen is in the top of the machine,
Easy for old people to read.

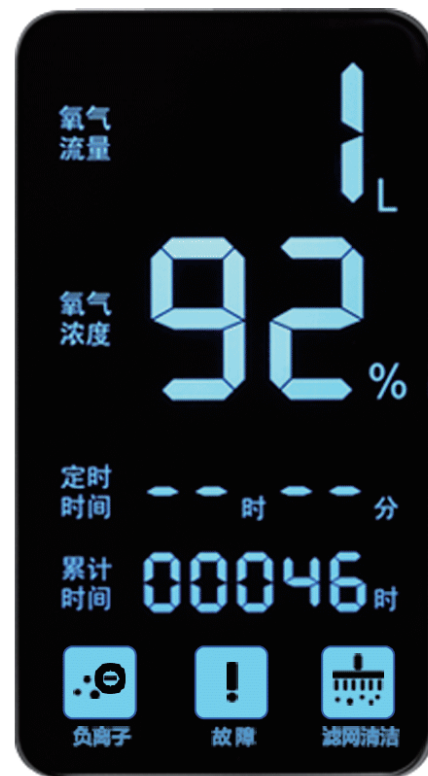

7.

Detachable filter

Convinient to replace or clean the filter.



8. Negative oxygen ions

Clean air, make negative oxygen ions


Make ozone and anion oxygen at the same.
Efficient to kill bacteria, air precipitation particles in the air.
Make the air more pure and fresh.



PSA

The filtered air will go into the compressor and it will compress the air into the molecular sieves, one of the molecular sieve produces oxygen with Pressure Swing Adsorption(PSA) and another starts low pressure desorption. The oxygen will then store in the oxygen tank before gives to patient.

Its environmental friendly and safe to use.

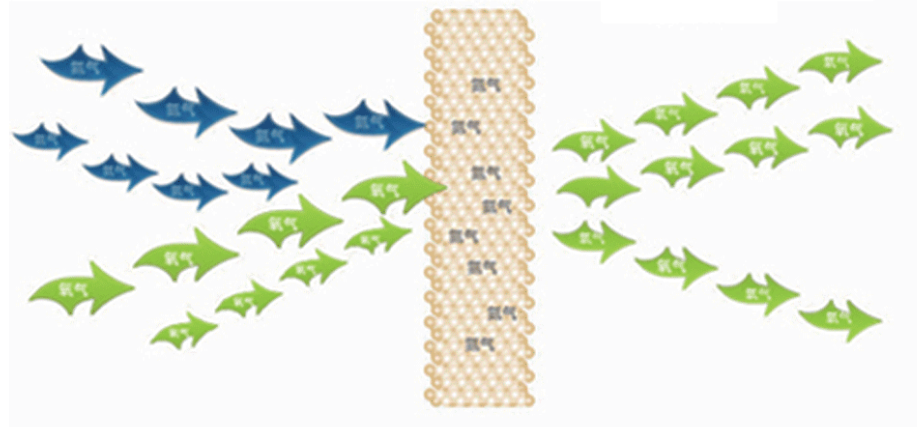

Life of compressor is more than 50000 hours



Yuwell homecare oxygen bar
Take it, it's yours!



Have you ever thought, in fact you are not sick just due to hypoxia caused by ever the
or to hypoxia caused by ever the hypoxia caused by ever the just due to hypoxia caused by ever thought, in fact you are not sick, it to hypoxia caused by ever this, in fact you are not sick just due to hypoxia caused by ever the
just due to hypoxia caused by ever thought, in fact you are not sick, it to hypoxia caused by ever the ever thought, in

YU300 Homecare oxygen concentrator

FRONTY VIEW >

YU300 HOME CARE OXYGEN

OXYGEN OUTLET

NO SMOKING

THE MAXIMUM
WATER LEVEL

开关

定时加

定时减

流量

负离子

POWER SWITCH


TIME INCREASING

TIME DECREASING

FLOWRATE

Negative oxygen ion



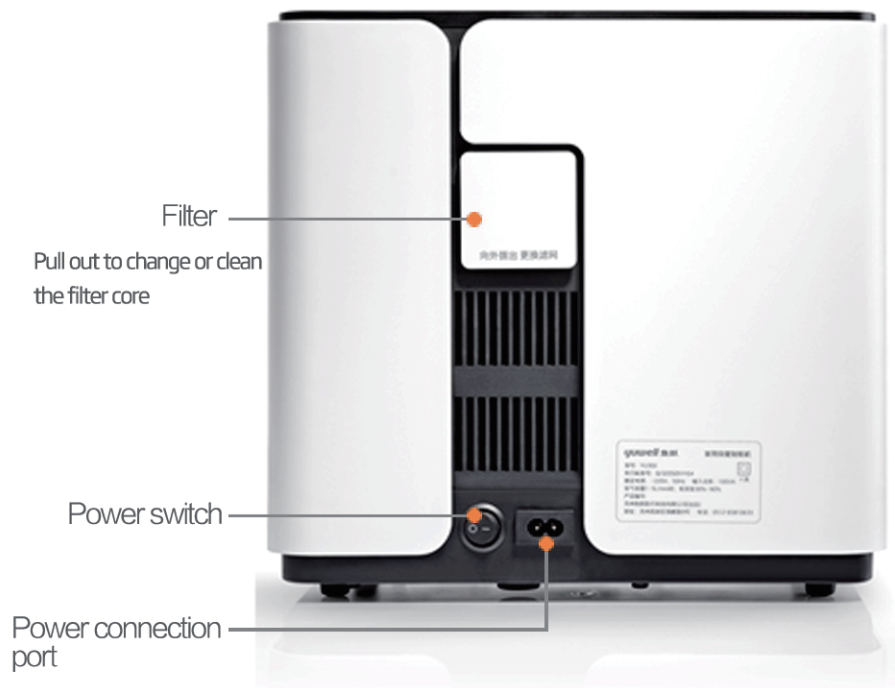

◀ SIDE VIEW





High
concentration
 $\geq 90\%$

Noise
 $\leq 40\text{dB(A)}$

Weight
only 8kg

Oxygen concentration	1~5L/min	Size	29.5 × 17.5 × 27 (cm)
Input power	100VA	Working system	Continuous
Electric classification	Class II	Power supply	AC220V ± 22V, 50HZ ± 1HZ
Safety system	Power failure alarm, compressor failure alarm, automatic shut off.		

Power cord



Nasal oxygen cannula



FREE

Filter core



FREE

Warranty card



User's manual



Product

